

croma

For creators of beauty.



saypha®
Hyaluronic acid
treatments

saypha®

Areas of application

for Croma HA fillers with lidocaine

Midface



Nasolabial folds



Lips




saypha® portfolio

saypha®
RICH

saypha®
FILLER

saypha®
VOLUME

saypha®
VOLUME PLUS

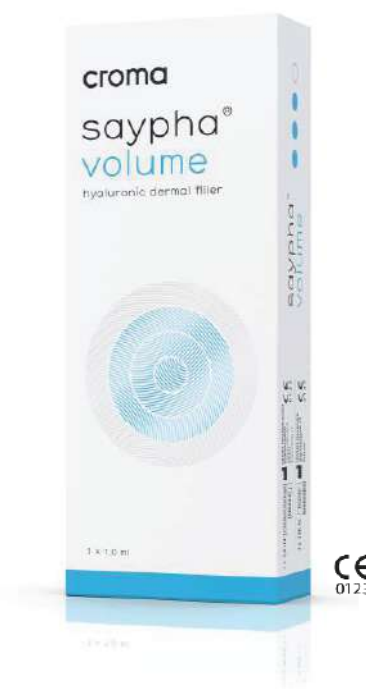

Product	HA concentration	Level of crosslinking	Mixing ratio crosslinked:non-crosslinked	Needles
saypha® RICH	18 mg/mL	□ □ □ □	-	2 x 30G ½" thin wall Terumo™ (CE 0197)
saypha® FILLER	23 mg/mL	■ ■ □ □	91:9	2 x 27G ½" thin wall Terumo™ (CE 0197)
saypha® VOLUME	23 mg/mL	■ ■ ■ □	95:5	2 x 27G ½" thin wall Terumo™ (CE 0197)
saypha® VOLUME PLUS	25 mg/mL	■ ■ ■ ■	95:5	2 x 27G ½" thin wall Terumo™ (CE 0197)

Hyaluronic acid + Lidocaine


saypha® FILLER Lidocaine

saypha® VOLUME Lidocaine

saypha® VOLUME PLUS

Croma syringe

Innovation for safety & comfort

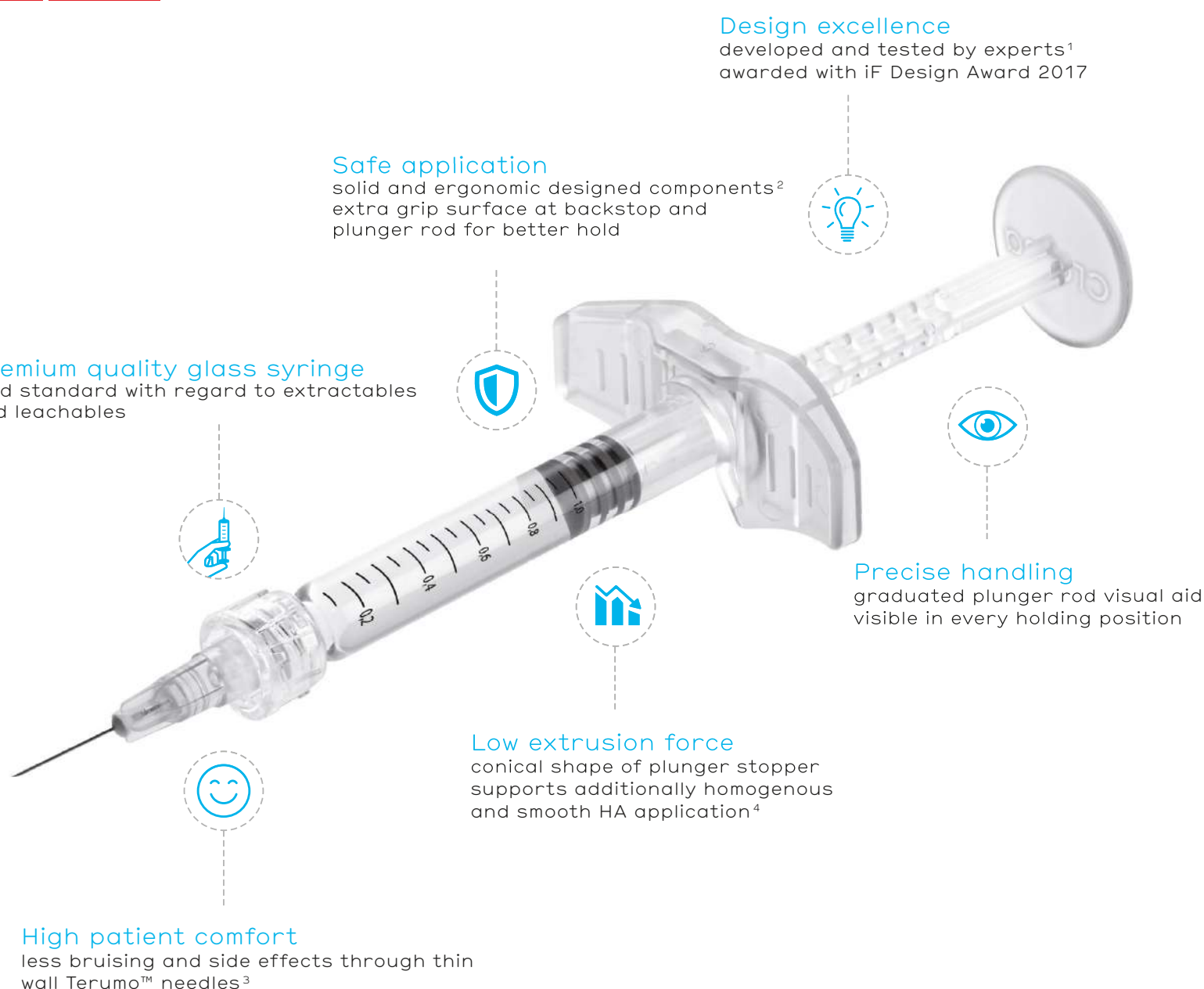

Extrusion force [N]

Product	EF [N] ⁵	Needle
saypha® FILLER*	10	27G 1/2" thin wall Terumo™ (CE 0197)
saypha® VOLUME*	12	
saypha® VOLUME PLUS**	15	

- A low extrusion force is beneficial for smooth and predictable injections.
- Croma HA fillers with lidocaine can be easily injected when using a 27G 1/2" thin wall needle.

Croma HA fillers

Crosslinked hyaluronic acid with 0.3% lidocaine

saypha® FILLER Lidocaine

Intended purpose	To create volume in order to correct wrinkles and folds in order to treat signs of ageing.
Indication	The device is indicated to correct moderate to severe nasolabial folds and to increase lip volume.
Injection area	Mid to deep dermis and submucosa
Ingredients	
Concentration HA	2.3% (23mg/mL) high molecular weight hyaluronic acid of non-animal origin
Crosslinking agent	BDDE (concentration ≤ 2 ppm)
Additional ingredients	Phosphate buffer, NaCl, 0.3% lidocaine hydrochloride
Packaging	
Packaging unit	1 box of 1 mL syringe
Needle	2×27G ½" thin wall Terumo™ needles (CE 0197)
Duration	
Est. duration in the skin	Anticipated up to 9 months ⁶

saypha® VOLUME Lidocaine

Intended purpose	To create volume in order to correct wrinkles and folds in order to treat signs of ageing.
Indication	The device is indicated to correct moderate to severe nasolabial folds.
Injection area	Deep dermis or subcutis
Ingredients	
Concentration HA	2.3% (23mg/mL) high molecular weight hyaluronic acid of non-animal origin
Crosslinking agent	BDDE (concentration ≤ 2 ppm)
Additional ingredients	Phosphate buffer, NaCl, 0.3% lidocaine hydrochloride
Packaging	
Packaging unit	1 box of 1 mL syringe
Needle	2×27G ½" thin wall Terumo™ needles (CE 0197)
Duration	
Est. duration in the skin	Anticipated at least 9 months ⁷

saypha® VOLUME PLUS

Intended purpose	The intended purpose of the device is to restore facial volume in order to correct moderate to severe midface volume deficit to treat signs of ageing.
Indication	The device is indicated to correct moderate to severe midface volume deficit in the zygomaticomalar region, anteromedial cheek region and submalar region.
Injection area	Deep subcutaneous and/or supraperiosteal
Ingredients	
Concentration HA	2.5% (25mg/mL) high molecular weight hyaluronic acid of non-animal origin
Crosslinking agent	BDDE (concentration ≤ 2 ppm)
Additional ingredients	Phosphate buffer, NaCl, 0.3% lidocaine hydrochloride
Packaging	
Packaging unit	1 box of 1 mL syringe
Needle	2×27G ½" thin wall Terumo™ needles (CE 0197)
Duration	
Est. duration in the skin	Approximately 12 months ⁸

Predictable results

thanks to highest quality and safety

The unique combination of high molecular weight hyaluronic acid of non-animal origin, a reliable production process in expert hands to crosslink the hyaluronic acid chains and a thorough purification process result in the exceptional safety profile of our Croma HA fillers.

Croma – the experts in hyaluronic acid

Croma has been manufacturing hyaluronic acid products for decades. As experts in HA with state-of-the-art facilities we pride ourselves on the high safety profile of our products. A glass syringe, low concentration of BDDE, stringent safety requirements and the overall performance of our Croma HA fillers are the reasons why doctors choose our aesthetic products.



croma
For creators of beauty.

Exclusive Distirbutor:

REDO
M A R K E T I N G

All right reserved by **Redo Marketing**